

PRODUCT INFORMATION
OTICON CONNECTLINE PHONE 2.0

Oticon | ConnectLine

*The **Oticon ConnectLine phone adapter** is an interface between a conventional landline phone installation and the Oticon Streamer, in effect turning a pair of Oticon hearing instruments into a wireless headset.*

As part of the ConnectLine system, the phone adapter works seamlessly and simultaneously with other ConnectLine solutions such as mobile phones, TV and the ConnectLine Microphone

A fully integrated system

With ConnectLine the user is in control of the system, not the other way around. Advanced features such as AutoPriority and AutoConnect ensure that the user can hear and answer a phone call while watching TV.

Landline home phone

Phone adapter 2.0 is installed in parallel to an existing landline home phone installation.

Once a user is in range of the phone adapter, the Streamer is automatically connected and ready to receive calls (AutoConnect).

The call can easily be answered, with a single touch of a button on the Streamer, at a distance of up to 30 meters from the phone adapter.

AutoPriority

If a user is listening to the TV using the ConnectLine TV solution (or other audio sources), the AutoPriority feature ensures that an incoming call is notified.

AutoResume

Once the call is ended the Streamer automatically resumes to its previous state. E.g. back to the sound from the TV.

For more information please refer to:
www.oticon.com/phone




PRODUCT OVERVIEW

Oticon | ConnectLine

Included

- RJ11 - RJ11 cable for connection to landline home phone
- Power supply
- Instructions for use

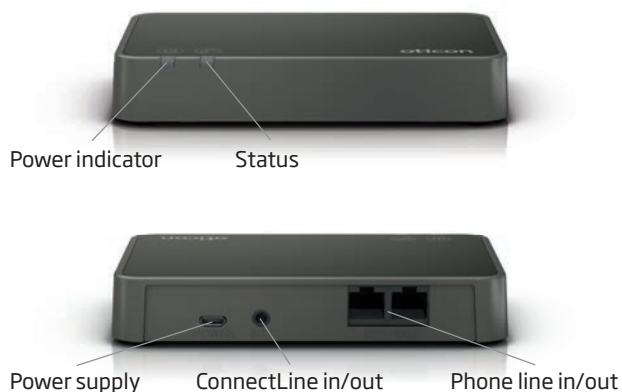

Technical data	
Compatibility	Streamer 1.4 and later
Dimensions	124 x 80 x 21 mm; weight 100 g.
ConnectLine system features	AutoConnect: Will automatically connect to the Streamer when in range. AutoPriority: Incoming calls will be prioritised over e.g. TV sound. Proximity pairing: The phone adapter will automatically pair to the Streamer if the Streamer is in pairing mode.
Range	Up to 30 meters with Streamer Pro (free line of sight)
Connections	Phone line in/out for interface to (analogue) landline home phone (PSTN) (RJ11/RJ11) ConnectLine in/out for direct interface to e.g. computer ¹
Power consumption	Max. 190 mA (using the included power supply)
Headset in	Max 6dBmV(RMS), 10 K Ω
Headset out	-15 – -10dBmV(RMS), 10 -100 Ω
Audio Bandwidth	License free 2,4 GHz (ISM band)
AutoConnect (range)	Yes (30 meters)
Audio Quality	Mono / 8 kHz sample frq
Frequency	80 Hz - 3.5 kHz (wireless)
Operating temperature	10 - 40 °C
Storage temperature	-20 - 70 °C
Emissions	Max. 39 dB μ V/m @ 3 meters. In compliance with all relevant standards

References:

1. A Computer interface cable is included with Streamer Pro. Also available as an accessory.